BellaBambi®

BellaBambi® multi

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BellaBambi® *multi*

Thank you for buying BellaBambi®
Congratulations on buying BellaBambi® multi. This is a multi-purpose and durable product for beauty, relaxation and sports. BellaBambi® consists of skin-friendly silicone and is produced exclusively in Germany. To have as much success and lasting pleasure with your product, we kindly ask you to read these instructions for use carefully.

Thank you!
Your team from BellaBambi®

SCOPE OF DELIVERY / DESCRIPTION / LEGEND TO SYMBOLS

Suction Pump

1 Turning knob
2 Tube connector
3 On/Off button

Application
BellaBambi® original

1 Tube connector

Application
BellaBambi® mini

1 Tube connector

Disinfectant spray

Flexible tube

Optional accessories

Application
BellaBambi® original

1 Tube connector

Application
BellaBambi® mini

1 Tube connector

Flexible Y-tube, to use two applications simultaneously

Optional accessories are available on www.bellabambi.de

WARNING This symbol warns against risks of injury or damage to your health

The device may only be used indoors

ATTENTION Indicates possible damages on the device/accessories

Insulated device (Protection class II)

INFORMATION Highlights important information

Not suitable for children under three years of age due to small partst

CORRECT USE AND FUNCTIONAL PRINCIPLE

BellaBambi® multi is to be used for pleasant and effective massages with negative pressure – on yourself or on other persons. Massages can be relaxing or stimulating, and are popular for treatment of stiff or tense muscles, pains and fatigue.

With BellaBambi® multi you can provide a negative pressure massage on individual parts of the human body. The device is not a medical product and is not indicated for medical purposes.


BellaBambi® is not to be used under any circumstances in the following cases:

- high-risk pregnancy
- on injured skin, skin dysfunctionality and skin inflammations, e.g. abrasion, ulcers, sunburn
- blood and lymphatic diseases
- skin cancer
- Do not apply BellaBambi® on:
 - Eyes and eyelids
 - Ears, lips and mouth
 - Nipples and sexual organs
 - Anus
- Known allergies against product materials and accessory components

In the following cases, the manufacturer recommends a comprehensive examination of the circumstances: pregnancy, complex regional pain syndrome (CRPS), cancer, tumors, disposition to hematomas and ulcers (e.g. caused by vascular diseases), irritable skin disorders such as psoriasis or urticaria (hives), with blood-thinning medication therapies, as well as any unclear medical conditions.

BellaBambi® causes a negative pressure on the skin, which uplifts the fascial tissue. The resulting space facilitates the flexibility of fascia and leads to tangible more smoothness. The innovative and unique peeling-rim generates higher shearing and a fresh and clear skin.

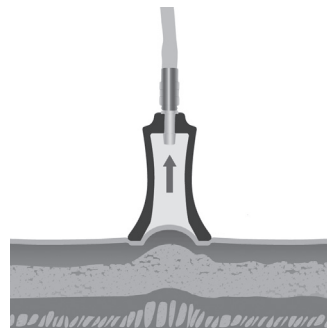Should your skin react sensitively, e.g. haematoma form quickly on the skin or if the skin becomes red, irritated or itchy, or general discomfort or pain is caused, then immediately discontinue with the treatment and consult your physician.

CORRECT USE AND FUNCTIONAL PRINCIPLE

i Please note our experiences with possible skin reactions when BellaBambi® is permanently left in place with medium suction power:

	pressure mark	redness	cupping mark*
sensitive skin	approx. 1-2 min	approx. 2-3 min	approx. 5 min
normal skin	approx. 2 min	approx. 4 min	approx. 7 min
robust skin	approx. 3 min	approx. 6 min	approx. 9 min

*cupping marks are not to be confused with haematomas or bruises which result from traumatic influences (e.g. punches, spills, etc.). Cupping marks indicate blood extravasation into the fascial tissue, which happens mostly without damaging vessels, and is never (!) painful. Furthermore, these marks disappear without becoming yellowish and greenish.

SAFETY PRECAUTIONS

! **Attention:** For your own safety, please adhere to the following safety precautions and make other users of the device aware of these. Keep this User Manual in a safe place.

- Before putting this device into operation, read the user manual carefully and completely and adhere to the sequence of the procedures.
- After usage suction pump and applications can be restored in the cardboard box. There is no packaging material that has to be disposed.
- This device must only be connected to alternating current. Check that the voltage indicated on the label of the casing corresponds to that of the mains voltage.
- Do not pull the power cord over sharp edges, nor pinch it or let it hang free. Keep the power cord well away from heat and oil.
- Ensure that the plug is close to the place where you wish to use this device.
- Avoid tripping over power cord and flexible tubes.
- The device may not be placed on the floor in the bathroom.
- Do not use near water.
- Do not immerse the device in water. Under no circumstances may water or other fluid penetrate the device!
- Do not cover the device with garments, towels or other materials to avoid overheating.
- Do not place the device on hot surfaces, such as stove plates, and do not use near an open flame.
- Switch of the device after each use. Do not pull on the power cord to remove the mains plug. Never plug in, unplug or otherwise handle the plug with wet hands.
- Do not use the device without supervision.
- Unplug the device before cleaning.
- Regularly check the device and the accessories for possible damage. The device may not be used if:
 - the casing of the pump or the power supply has been damaged
 - it has been dropped onto a hard surfaceIn both instances, qualified technical personnel or our service department should inspect the device. Never try to repair the device yourself.
- Children after the age of eight (8) and persons with restricted physical, sensory or mental capabilities and/or by persons lacking the required understanding of the equipment and how it is used, unless they are supervised by a person who is responsible for their safety or they are given instructions on how to use the device. Children should be supervised to ensure that they do not use the device as a toy. Children are not allowed to clean the device without supervision.
- Do not exceed dosages and time limits as recommended in this user manual.
- Only use accessories that have been recommended by the manufacturer.

1. Before each usage clean all applications as described in section „Cleaning“ on ► page 23
2. Select the appropriate BellaBambi® application according to the size of the body part and connect it to the suction pump with the flexible tube.
 - BellaBambi® mini is suitable for small and filigree parts of the body, e.g. face, hands and feet.
 - BellaBambi® original is most suitable on larger surfaces, e.g. back, legs and stomach.
3. Plug in the device and switch on.
4. Turn the knob to the maximum right and place BellaBambi® on your skin. As soon as it has adhered, reduce the suction power until it is recognizable and comfortable if necessary
5. Please follow the recommended steps on your desired application:
 - Cellulite ► page 19
 - Peeling and skin care ► page 19
 - Fascial fitness ► page 20
 - Relaxation ► page 20
 - Scar treatment ► page 21
 - Reflexology ► page 21
 - Massages ► page 22
6. Clean thoroughly after each use, as it is described in section „Cleaning“ ► page 23
7. Switch off and unplug the device from the mains socket.

Please note that the suction pump is equipped with an internal electronic timer, which automatically counts down and switches off the pump after 15 minutes. The countdown always starts anew when the pump is switched on.

APPLICATIONS


Always use BellaBambi® on skin- and muscle tissue. Never apply it on the spinal cord or on protruding bones. Glands and lymph nodes may not be treated.

Always adopt the suction power to your personal well-being.

CELLULITE


Negative pressure massage of connective tissue is vitalising and refreshing too. BellaBambi® can be combined with cellulite creme and body lotions to improve your skin texture.

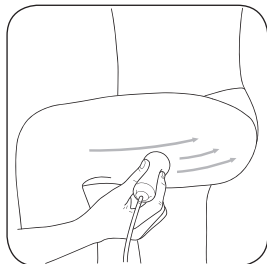Place BellaBambi® above your knee and drag it to your hip. Relocate BellaBambi® until your entire thigh is covered.

Repeat 3 times altogether.

Tip: Do sports before or after using BellaBambi®, especially swimming or cycling. Workouts for legs and buttocks suit perfectly too.

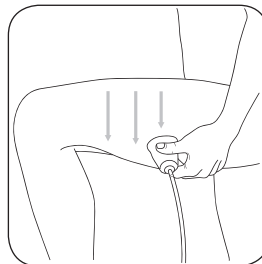Place BellaBambi® approx. 5 cm above your knee and drag it around your thigh without stopping. Relocate BellaBambi® above and repeat the treatment until your entire thigh is covered.

Repeat 3 times altogether.

Tip: Make sure to drink enough water or tea and eat fresh vegetables. Avoid sugar and alcohol where possible.

PEELING AND SKIN CARE

BellaBambi® revives and softens your skin and does not damage the environment through microplasticity. The unique peeling-rim allows increased shearing and leaves the skin looking smooth and fresh.

Place BellaBambi® below your hair-line in the middle of your forehead and drag it outward. Relocate BellaBambi® 1 cm below and repeat until you reach your eyebrow. Move to your cheek and drag BellaBambi® from the side of your nose to your ear.

Important: Keep sweeping BellaBambi® to minimize the risk of cupping marks.

Place BellaBambi® in the middle of your chin under your bottom lip and drag it over your jaw to your earlobe. Repeat until you've finished your chin. Finally place BellaBambi® below your jawbone and drag it from your upper neck to your collarbone.

Important: Please do not treat your larynx.

BellaBambi® causes a negative pressure on the skin, which uplifts the fascial tissue. The resulting space facilitates the flexibility of fascia and leads to tangible more smoothness.

Place BellaBambi® below your knee, drag it to your wrist and back again. Relocate BellaBambi® until your entire calf is covered. Repeat 3 times.

Variation:

Try to stretch and bend your foot while dragging BellaBambi®.

Tip: On the sides of your knee, you can also drag BellaBambi® until your thigh.


Lumbar fascial tissue is the biggest fascial complex in the human body and often causes problems.

To preserve its smoothness place BellaBambi® aside your spine and drag it as the arrows show.

You can apply fascial fitness for 10 minutes a day.

RELAXATION

Relaxed well-being is a welcome side effect, which arises already after a short period of time. The negative pressure not only lifts your skin, but also your mood.

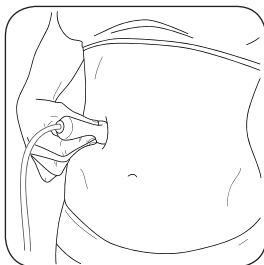
Place BellaBambi® with low suction power for 20 seconds above your solar plexus. It is located on half of the distance from your breastbone to your belly button.

Important:

Focus on that point and breathe easy, deeply and slowly.

Tip: Another person could help you to replace BellaBambi®, so you can concentrate on focusing and relaxing.

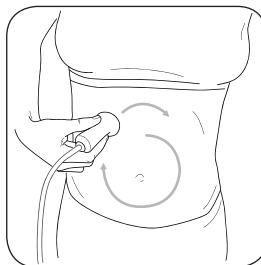Place BellaBambi® below your costal arch and drag it slowly in narrowing circles around your belly button.


Afterwards do it the reverse way until you have reached the initial point. This exercise is best done lying down.

Important:

Wind down and breathe easy.

SCAR TREATMENT

Stretching tridimensional can release fascial adhesion. Treatment with BellaBambi® can improve your quality of life through greater flexibility.

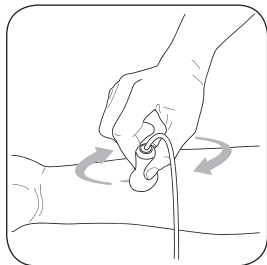
Warning:

Please consult your physician before using BellaBambi® on scars! Be careful!

Important:

The scar has to be clean, dry and closed.

Place BellaBambi® on the scar, grab it on the top and circuit. Alternatively, push BellaBambi® in different directions for 20 seconds each.

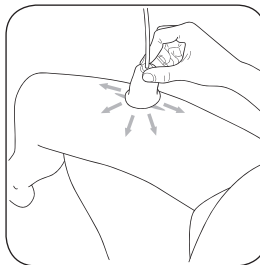Place BellaBambi® and drag it in a star pattern over your scar. Please fix the initial point with your other hand, when dragging BellaBambi®. Start with low suction power and don't use oil or other lubricants.


Warning:

Please consult your physician before using BellaBambi® on scars! Be careful!

Important: The scar has to be clean, dry and closed.


REFLEXOLOGY

You can precisely apply BellaBambi® on your sole of foot. The negative pressure releases the pressure, you literally are „under“ the whole day.

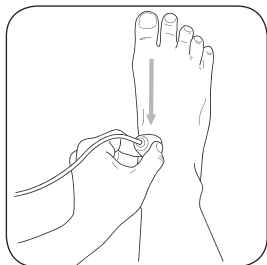Place BellaBambi® below your toes and drag it to your ankle. You can also set it between your toes' extensions. Repeat on the bottom of your foot and drag BellaBambi to your heel.


You can use massage-oil or other lubricants to facilitate gliding.

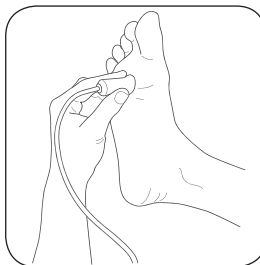You can work lengthwise and crosswise your sole of foot to fulfill a stimulating reflexology with BellaBambi® mini, or also with BellaBambi® original, in case you have big feet.

Treating for 5 minutes should be sufficient on each foot.


BellaBambi® prepares your body with a pleasant and warming massage. If you drag BellaBambi® on your skin it is even self-massaging, and it preserves your masseur's joints.

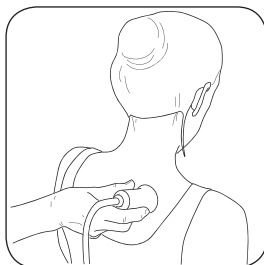Place BellaBambi® below your hairline and drag it downwards as far as possible. Afterwards relocate BellaBambi® to the side and repeat.

A manual follow-up massage is perceived as particularly beneficial – even the masseur appreciates the preparatory work of BellaBambi®.

Attention: Do not use on the spine.

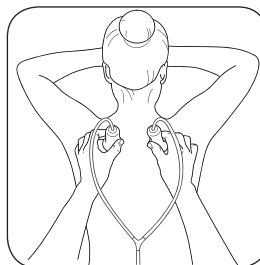If used on another person's back, drag BellaBambi® crescent-shaped outwards.

Speed up or slow down according to your individual preferences.

Attention:
Do not use on the spine and always listen to your body's signals.

i Flexible Y-tube and other BellaBambi® applications are optional accessories.

Optional accessories

Application
BellaBambi® original

Application
BellaBambi® mini

Flexible Y-tube, to use two applications simultaneously

Optional accessories are available on www.bellabambi.de

CLEANING

Clean the device only with a moist (never dripping wet!) cloth. If necessary, use a few drops of dishwashing detergent.

1. Hand wash BellaBambi® applications in warm water using mild soap or diluted dishwashing detergent, or place it in dishwasher.
2. Dry BellaBambi® application thoroughly.
3. Disinfect BellaBambi® by applying a disinfectant spray, suitable for silicone. The disinfectant should not contain phenol derivatives or alkylamines (glucoprotamines). We recommend: Desomed® Rapid AF, which is available on www.bellabambi.de.
4. Leave disinfectant on for a while, according to the instructions of the manufacturer.
5. Then rinse again with water to remove any residues.
6. Allow BellaBambi® to dry thoroughly.

IN THE EVENT OF PROBLEMS

PROBLEM	CAUSE	REMEDY
BellaBambi® does not adhere to the skin	Suction pump is not connected to the power supply or switched off.	Connect the plug to the mains supply and switch on the appliance.
	Applications are not connected correctly.	Control the connecting parts of the suction pump and the application.
	Skin is very hairy.	Use some lubricant, eventually shaving could be necessary.
	The opening area of BellaBambi® does not cover completely.	Use a smaller BellaBambi® or squeeze it a little bit.
	Timer has expired.	Use the On/Off button.
Application is unpleasing	Suction power is set too high.	Reduce suction power until it is noticeable but comfortable.
	No lubricating skin care is used.	Use lubricating skin care, if possible.
	Not enough lubricating skin care is used.	Use more lubricating skin care, or another product, which does stronger support gliding.

You can contact BellaBambi® if you have any questions. By phone, e-mail or facebook:
+49 (9073) 921 046 1 | info@bellabambi.de | [facebook.com/bellabambi.de](https://www.facebook.com/bellabambi.de)

TECHNICAL DATA AND AMBIENT CONDITIONS

The permitted ambient conditions define the limits within which BellaBambi® multi can be used without danger of damage

Permitted ambient conditions:

- +10°C to +40°C (in use)
- -18°C to +60°C (storage)
- 30% to 75% relative humidity

Power consumption: 230 Vac, 50 mA, 50/60 Hz

	Adjustment	Negative pressure at device outlet
	Intensity min (turning knob maximum left):	-140 mbar
	Intensity max (turning knob maximum right):	-280 mbar

MAINTENANCE AND REPAIR

Any repairs on the suction pump must be carried out by the manufacturer only. Never open the device yourself. Please do not manipulate the device in any way. Otherwise, the warranty will become null and void. Should your device be defective, please contact the following address:

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DISPOSAL

Your device contains various materials, which can be recycled. To dispose of the device, take it to a local municipal waste collection point. Do not dispose of the device with your household waste. Please enquire at your municipal recycling facilities of possibilities for throwing away old or defective electronic equipment. If this is not possible, you can send the device back to the manufacturer for proper disposal. In this case, the sender will pay postage or shipping costs.

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